

**MIDMARK®**

# Installation Instructions

## 9A278002 Wrist Rest Kit

### KIT COMPONENTS

| Part Number  | Description              | Qty. |
|--------------|--------------------------|------|
| 029-2424-XXX | Wrist Rest Assembly..... | 2    |
| 003-1072-00  | I / O Instructions.....  | 1    |

### REQUIRED TOOLS

Thin blade screwdriver  
3/8" Wrench  
1/8" diam. Punch

### INSTALLATION


### INTRODUCTION

To provide support for wrists during examination and general procedures.

#### NOTE

9A278002 Wrist Rest Kit can only be used on Models 419-003, 004, 005, and 006 with the threaded mounting holes in the headrest housing.

For units without mounting holes use Wrist Rest Kit (002-0647-00).

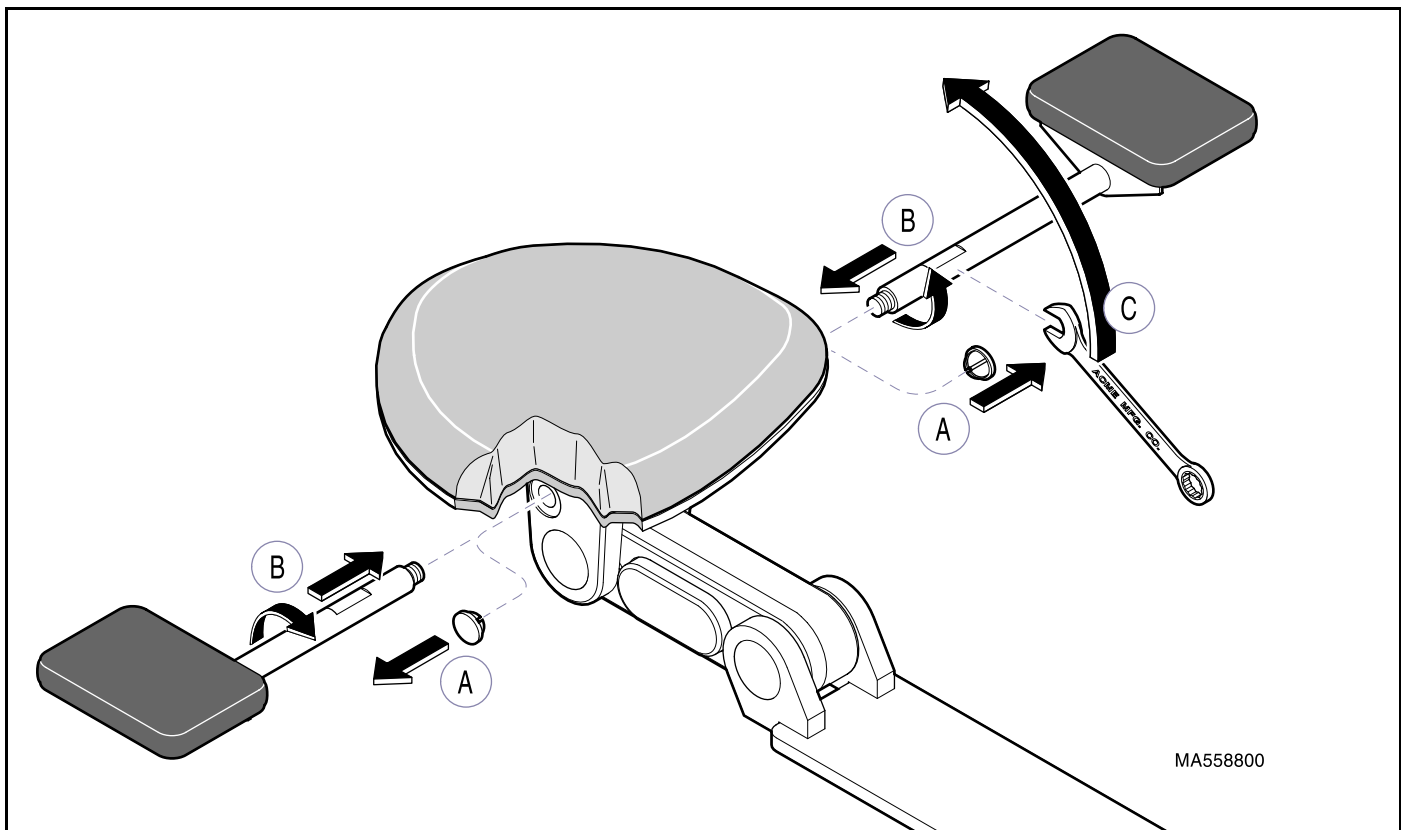

- Remove hole plugs from head rest.
- Hand thread each wrist rest assembly into mounting holes.
- Tighten wrist rest assembly using a 3/8" wrench.

**MIDMARK®**

Midmark Corporation, Versailles, Ohio 45380 U.S.A.  
(937) 526-3662 FAX (937) 526-5542

003-1072-00 Rev. A1 (2/00)